


CHECK&TRACK VERSION 2.0

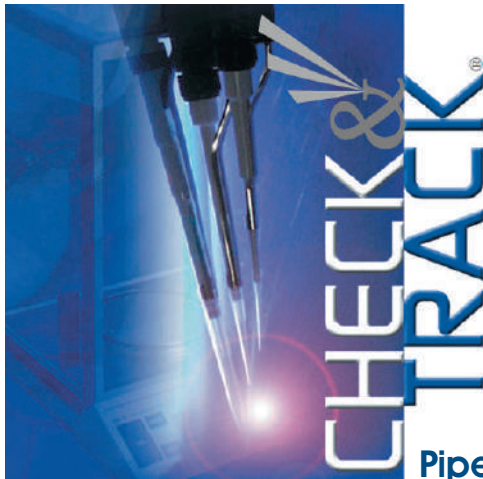

Pipette Calibration Software Version 2.0

Pipettes, the main volumetric instruments in the lab, are the core of your laboratory quality management. Therefore just-in-time acquisition of your pipettes' precision and accuracy is a main issue in a GLP or GMP facility. Check&Track® is a safe and convenient solution managing certification and record keeping of your pipettes according to the different QA standards.

Meets ISO 9000, ISO 8655, DIN 12650, NCCLS, GMP and GLP requirements:

Main Features of Check&Track®:

Check&Track® is suitable for:

System requirements:

- Flexible report layout including specific company information
- Complete and printable audit trail
- Specific administrator and operator functions
- Login and password administration
- Guided pipette checks
- Electronic collection of balance weighing values from any balance connected via serial port (RS232).
- Electronic collection of environment parameters from an environment sensor connected via USB (temperature, pressure and humidity).
- Automatic conversion of balance weighing values to volumes, taking into account temperature, pressure, humidity, altitude and balance weight density.
- Automatic calculation of statistical significance, based on the specified number of measurements, accuracy and precision.
- Complete pipette check history, including raw data.
- Accuracy and precision history (charts) for every pipette.
- Email reminders per operator, listing all pipettes due for checking within a configurable number of days.
- Printable audit trails.
- Single and network installation supported.
- IQ/OQ documents available on request.
- Fixed volume pipettes
- Variable volume pipettes
- Multi-channel pipettes
- Stepping pipettes
- Volumetric glassware
- Windows XP Service Pack 3 or later
- Screen resolution: 1280 x 1024

Download Check&Track® trial version from www.check-and-track.com

Do you have questions about Check&Track®? Ask our Application Specialist.